


CHROM-LINK®

MONITORING DATA ACQUISITION SYSTEM

CHROM-LINK® is a complete instrument control and monitoring data acquisition system dedicated to a single MINICAMS® monitor. Operating under Windows®, CHROM-LINK collects all MINICAMS data, including concentration and status reports, chromatograms, calibrations, and operating conditions into a relational database. This allows for a wide variety of reports to produce a complete audit trail of chemical-monitoring activities and instrument performance. CHROM-LINK performs reliable backup and archiving of air-monitoring and chromatographic data automatically, either on a regularly scheduled basis or on demand.

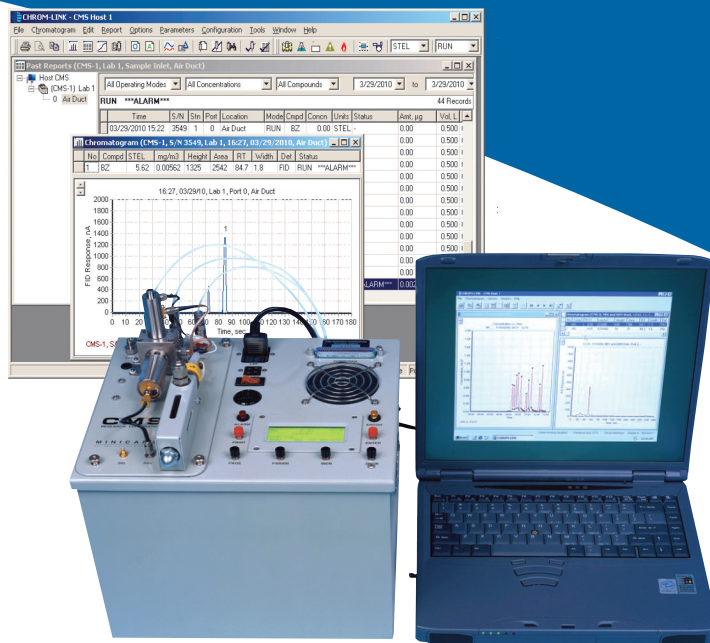


MINICAMS® plus CHROM-LINK®

FEATURES

- Chromatogram collection
- Reporting and auditing capabilities
- Multi-level calibration and curve-fitting options
- Remote instrument control
- Library of established methods
- Monitoring setup and configuration
- Software package included on laptop platform
- Package includes Multi-Level Calibration, Analytical Methods Library, simulated LCD screen, and other MINICAMS control functions
- Compound history view allows tracking of MINICAMS performance over time, including detector sensitivity, calibration frequency, and analytical standard consistency
- Logged data is easily retrievable and can be copied to other applications
- Off-line review of archived data with optional CHROM-LINK View software
- Standard navigation tools allow easy transition from concentration report screen to corresponding chromatogram, operating parameters, or other data views
- Filtering features allow preparation of concise event reports


Concentration Trend Graph

- User-selected ranges
- Menu-selected options include:
 - Display data values
 - Add annotation to a movable annotation field
 - Modify scales
 - Select compounds to be displayed

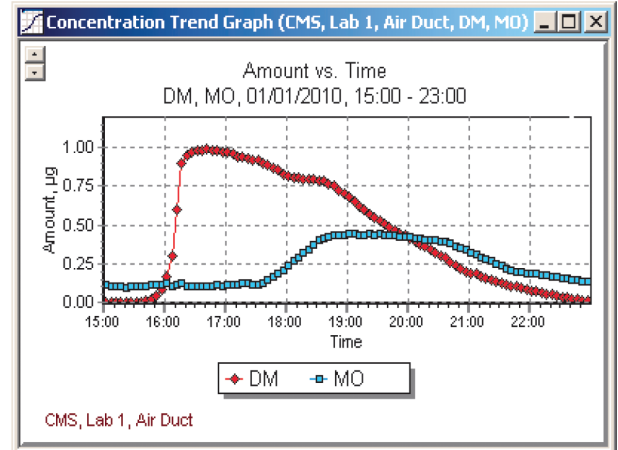


Calibration

- First-order, second-order, third-order, and power-fit calibration curves available
- Past calibrations retained in CHROM-LINK[®] database
- Calibration reports available for all calibrations

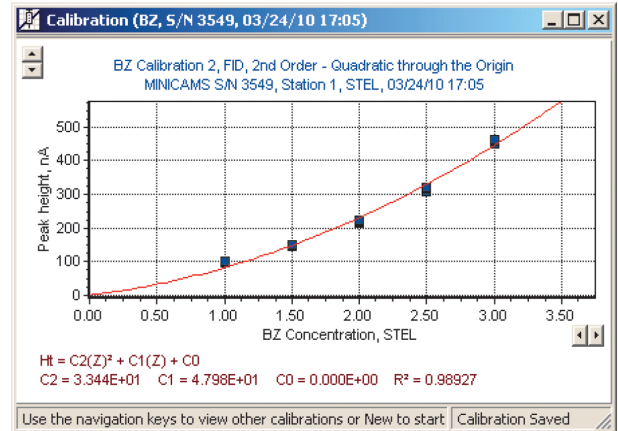

Chromatogram

- Chromatograms include concentration report information and can be navigated in chronological order

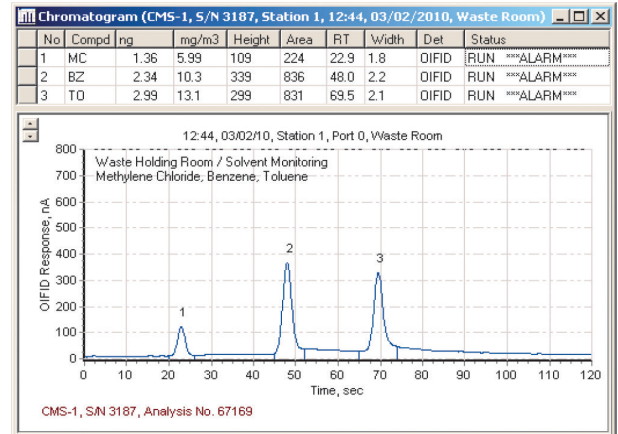

Contact OI/CMS for additional information:



2148 Pelham Parkway
Building 400
Pelham, Alabama 35124

(205) 733-6900
FAX (205) 733-6919

www.oico.com
E-mail: CMSSales@Xyleminc.com

MINICAMS and CHROM-LINK are registered trademarks of OI Analytical/CMS Field Products, Pelham, Alabama, USA.

Windows is a registered trademark of Microsoft Corp.

Publication 44310917